

NEW

Nichipet[®] *Air*

Autoclavable
MADE IN JAPAN
NICHIRYO

Autoclavable
MADE IN JAPAN
NICHIRYO

Features of *Nichipet[®] Air*

121°C
20min
Autoclavable



MADE IN
JAPAN

The lighter weight and reduced plunger pressure of the Nichipet Air helps mitigate fatigue and stress from repetitive pipetting, while still delivering world-class Accuracy and Precision.

1 Easy calibration

Easy-Calibration function provides simple lab calibration. (calibration jig needed)

2 Pipetting

Up to 70% reduction in plunger pressure compared with conventional models.

3 Finger rest



Ergonomically shaped to fit your hand.

4 Main body

As durable as previous Nichipet models but 56% lighter - only 70g!
Fully autoclavable.

5 Ejector




Removable ejector allows unit to be used in thin test-tubes.

Cat. No	Volume range (μL)	Increments (μL)	Color
DS-2-NAR	0.1~2	0.002	Grey
DS-10-NAR	0.5~10	0.01	Red
DS-20-NAR	2~20	0.02	Yellow
DS-100-NAR	10~100	0.1	Orange
DS-200-NAR	20~200	0.2	Light Blue
DS-1000-NAR	100~1000	1	Dark Blue
DS-5000-NAR	500~5000	5	Green
DS-10ML-NAR	1000~10000	10	Light Green



6 Smooth volume setting

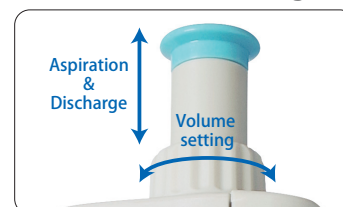

Volume setting rotator is independent from the push button, preventing volume shift during pipetting

7 Volume indicator window

The volume can be checked during pipetting.

8 Pipette volume

*Color key able at lower left.




Each volume has a unique color!



Three ways to easily recognize your pipette's volume.

9 Hyper blower system

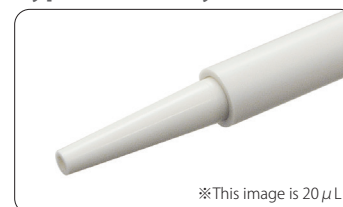

※This image is 20 μL

In order to completely discharge the set sample amount, the 2 μL, 10 μL and 20 μL units have a high-performance blow-out function.

Distributor

dot scientific inc

NICHIRYO CO.,LTD.

TOKYO MAIN OFFICE

Sakura Building 3F 1-10-1 Kandanishiki-cho, Chiyoda-ku,
Tokyo, 101-0054, Japan

TEL. +81-3-6273-7652

✉ info@nichiryo.co.jp www.nichiryo.co.jp/en/