




The Doble QWick Series utilizes an absorbent wicking material that charges with liquid nitrogen in fewer than two hours, providing the capacity for same-day vapor shipping. A unique absorbent layer in the base of the storage tanks enables them to be charged with nitrogen and employed as dry shippers with hold times of up to 30 days. Once at the final destination, the tanks can be filled with liquid and used for long term storage, therefore avoiding the need for return shipments. Equipped with the same capacity and dimensions as the comparable Doble Series.

## Features Include:

- CHARGES IN LESS THAN TWO HOURS
- Liquid and vapor storage options
- Low liquid nitrogen consumption
- Convenient lightweight package

## Accessories:



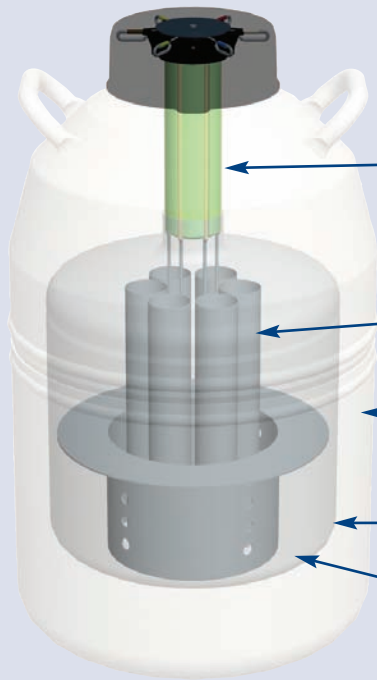

- A** Durable composite lid design
- B** High strength neck tube reduces liquid nitrogen loss
- C** Canisters
- D** Strong, lightweight aluminum construction
- E** Advanced chemical vacuum retention system
- F** Hydrophobic liquid nitrogen absorbent system

	Doble QWick 10/660	Doble QWick 20/660	Doble QWick 17/2400	Doble QWick 20/2100	Doble QWick 14/4500	Doble QWick 14/3500
<b>Comparable To</b>	Doble 11	Doble 20	Doble28	Doble 34	Doble 47	Doble 47-10
<b>Maxium Storage Capacity</b>						
<b>Number of Canisters</b>	6	6	6	6	6	10
<b>No. of 1/2 cc Straws 10/cane</b>	660	660	2400	2100	4500	3500
<b>Vial Capacity</b>	210	210	810	630	1320	1050
<b>Performance</b>						
<b>Liquid Nitrogen Capacity L</b>	10.0	18.5	28.0	32.0	46.0	46.0
<b>Vapor Capacity L</b>	3.1	2.8	8.4	6.7	8.4	9.6
<b>Static Evaporation Rate L/day</b>	0.17	0.10	0.35	0.20	0.40	0.40
<b>Normal Working Duration days for liquid</b>	37	116	50	100	74	72
<b>Normal Working Duration days for vapor</b>	10	20	17	20	14	14
<b>Unit Dimensions</b>						
<b>Neck Opening in (mm)</b>	2.18 (55)	2.18 (55)	3.81 (97)	3.50 (89)	5.00 (127)	5.00 (127)
<b>Overall Height in (mm)</b>	21.6 (549)	25.7 (652)	22.0 (559)	26.6 (676)	26.5 (673)	26.5 (673)
<b>Outside Diameter in (mm)</b>	10.2 (260)	14.5 (368)	18.2 (462)	18.2 (462)	20.0 (508)	20.0 (508)
<b>Canister Height in (mm)</b>	11 (279)	11 (279)	11 (279)	11 (279)	11 (279)	11 (279)
<b>Canister Diameter in (mm)</b>	1.65 (41)	1.65 (42)	3.09 (78)	2.81 (72)	4.00 (101)	2.81 (72)
<b>Weight Empty lb. (kg)</b>	17 (7.7)	25 (12.0)	34 (16.0)	39 (18.0)	47 (21.3)	48 (21.8)
<b>Weight Charged vapor lb. (kg)</b>	22.5 (10.2)	30.2 (14.0)	47.0 (21.3)	50.0 (22.7)	63.0 (28.6)	63.0 (28.6)
<b>Weight Full liquid lb. (kg)</b>	36 (16.0)	60 (27.0)	84 (38.0)	95 (43.0)	125 (57.0)	120 (54.5)

**THREE Year Vacuum Warranty**

**Chart BioMedical**  
 2200 Airport Industrial Dr., Ste. 500  
 Ball Ground, GA 30107  
 Ph 770-721-7759 • Toll Free 1-800-482-2473  
 Fax 770-721-7758  
 customerservice.usa@chartindustries.com



www.chartbiomed.com

**Chart BioMedical Ltd.**  
 Ph +44(0) 1344 403100 • Fax +44(0) 1344 429224  
 customerservice.europe@chartindustries.com  
**Chart Japan Co., Ltd.**  
 Ph + (03) 5776 2670 • Fax + (03) 5776 2676  
 customerservice.japan@chartindustries.com  
**Chart Australia Pty Ltd.**  
 Ph + 61 2 974 94 333 • Fax +61 2 974 94 666  
 customerservice.australia@chartindustries.com

Copyright © 2010 Chart Industries Ref ML130002 Rev A